


## Organic solvents

Data No. GA086-0000

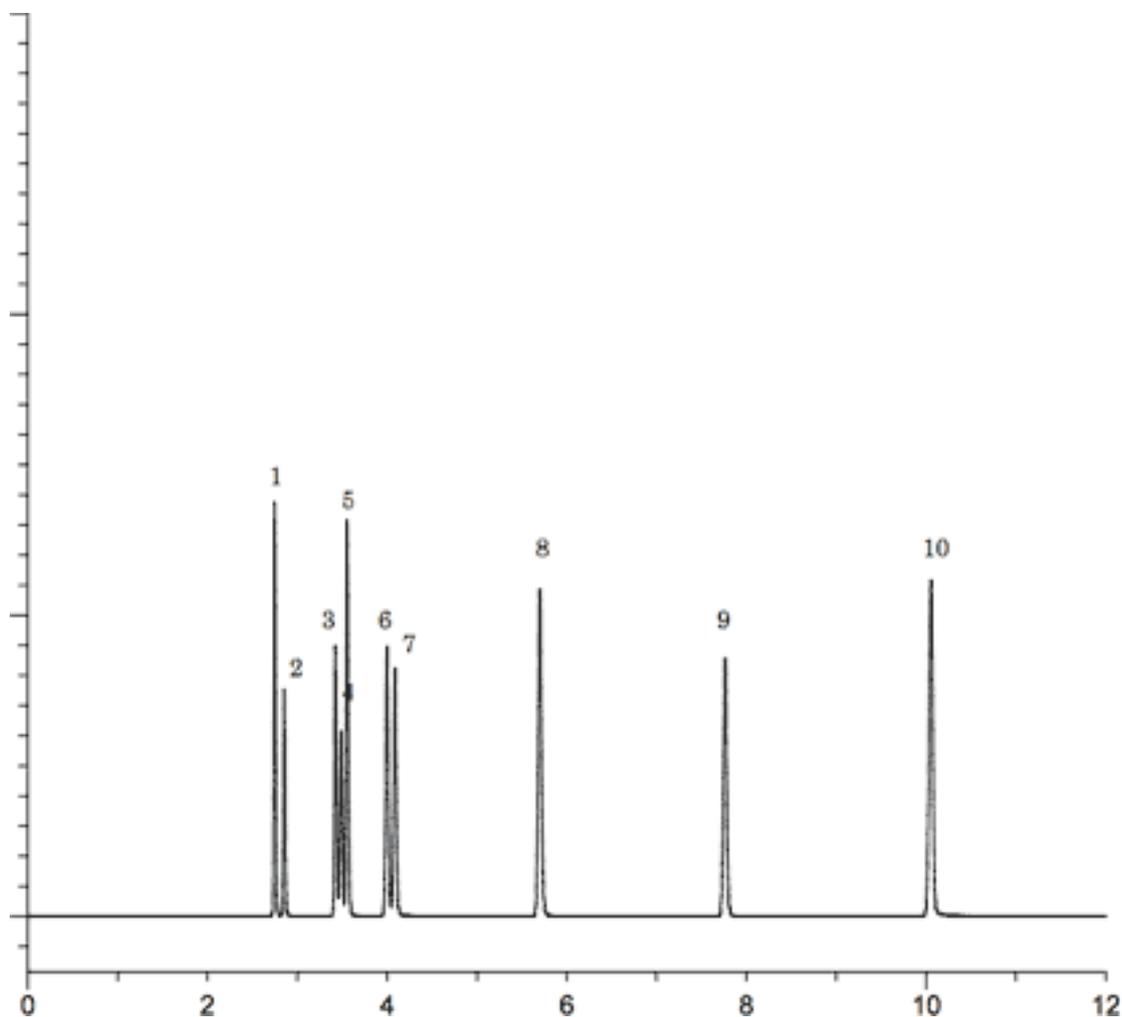

System : GC/FID  
Column : InertCap Pure-WAX  
0.25mm I.D. x 30m df=0.25µm  
Column Temp. : 40°C (5min) → 5°C/min → 100°C  
Carrier Gas : He 100kPa  
Injection : Split 1 : 80 0.1µl  
250°C  
Detection : FID Range 10°  
260°C  
Sample : Organic solvents

Solutes : 1. Acetone  
2. Methyl acetate  
3. Ethyl acetate  
4. Methanol  
5. Methyl ethyl ketone  
6. iso-Propanol  
7. Ethanol  
8. Methyl isobutyl ketone  
9. n-Butyl acetate  
10. n-Butanol