


## Purification of *Coptidis rhizoma*

Data No. LA145-0000

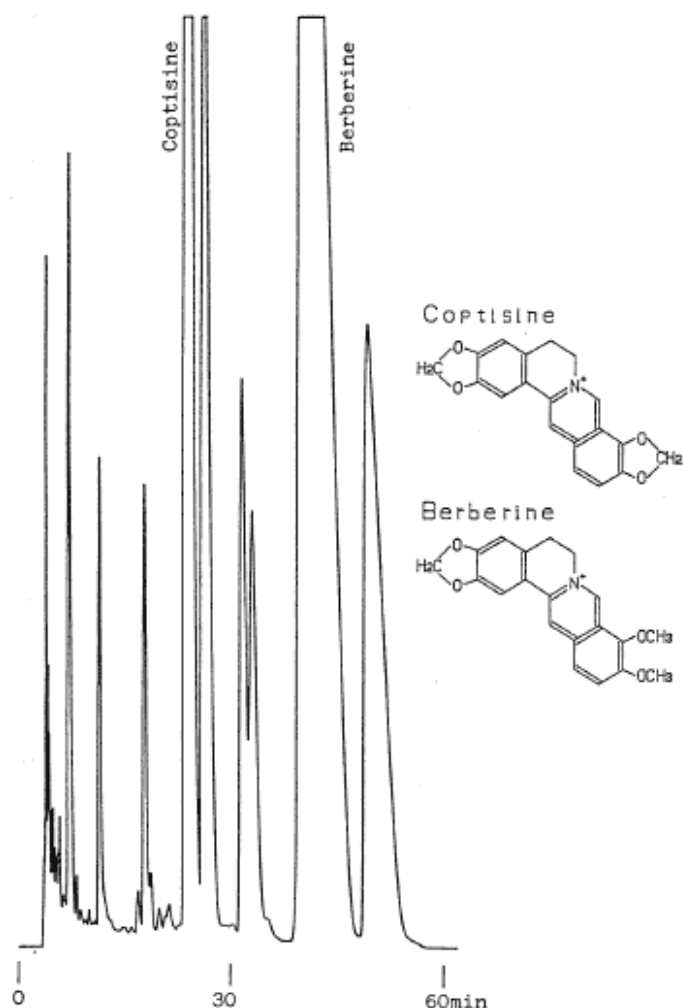

Column : Inertsil PREP-ODS  
10µm  
250 × 30mmI.D.

Solutes : 1. Coptisine  
2. Berberine

Column Cat. No. : 5020-31133

Eluent : CH<sub>3</sub>OH/0.1% CF<sub>3</sub>COOH=30/70

Flow rate : 30 mL/min.

Col. Temp. : 40°C

Detection : UV 254 nm

Injection Vol. : 2.5mL

Sample : Coptidis Rhizoma(50mg/mL)