


Analysis of Steroid Hormones

Data No. LA253-0000

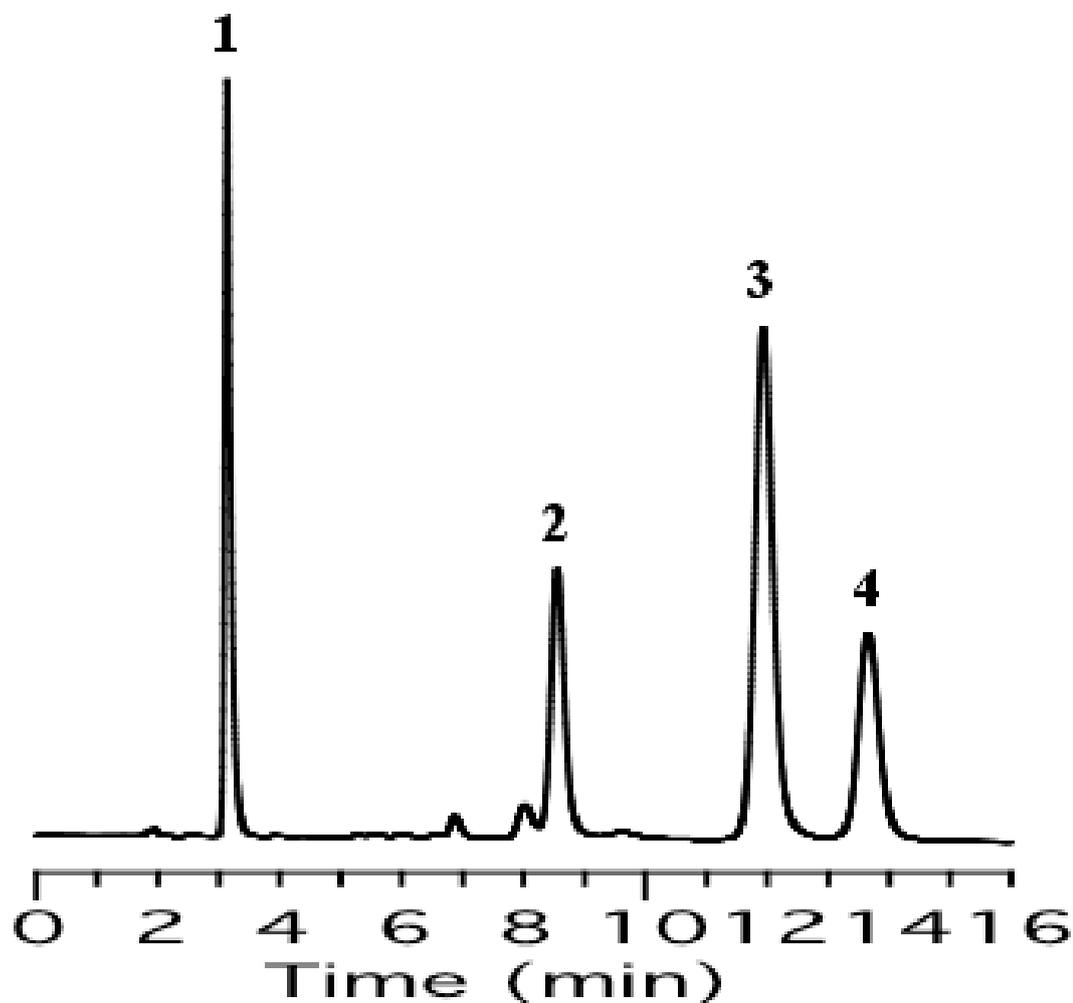

| | | | |
|-----------------|---|---------|--|
| Column | : Inertsil CN-3 5 µm 150 x 4.6 mm I.D. | Solutes | : 1. Progesterone 2. Corticosterone 3. Prednisone 4. Prednisolone |
| Column Cat. No. | : 5020-01940 | | |
| Eluent | : A) n-Hexane B) Ethanol A/B = 90/10, v/v | | |
| Flow rate | : 1.0 mL/min. | | |
| Col. Temp. | : 40°C | | |
| Detection | : UV 254 nm | | |
| Injection Vol. | : | | |
| Sample | : | | |