


# InertSearch™ for GC

InertCap® Applications

## Pesticides in golf courses

Data No. GA008-0000

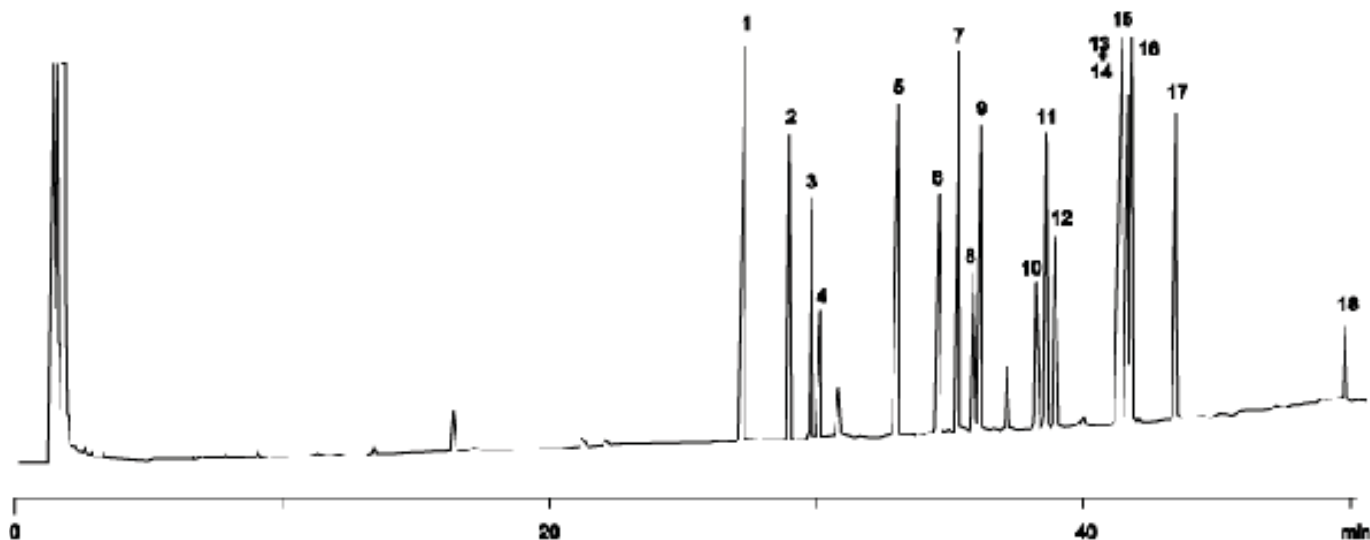

System : GC/FID  
Column : InertCap 5  
0.25mmI.D.×30m df=0.4µm  
Column Temp. : 100° C→3° C/min→250° C  
Carrier Gas : He 30cm/sec  
Injection : Split 1 : 70  
Detection : FID Range 1×8  
Sample : Pesticides

Solutes : 1. Simazine(CAT)  
2. Propyzamide  
3. Diazinon  
4. Chlorothalonil(TPN)  
5. Tolclophosmethyl  
6. Fenitrothion(MEP)  
7. Thiobencarb  
8. Fenthion(MPP)  
9. Chlorpyrifos  
10. Pendimethalin  
11. Captan  
12. Isofenphos  
13. Napropamide  
14. Butamifos  
15. Flutoluanil  
16. Isoprothiolane  
17. Isoxathion  
18. Iprodione