


# InertSearch™ for GC

InertCap® Applications

## 2-Pyrrolidone

Data No. GA 269-0904

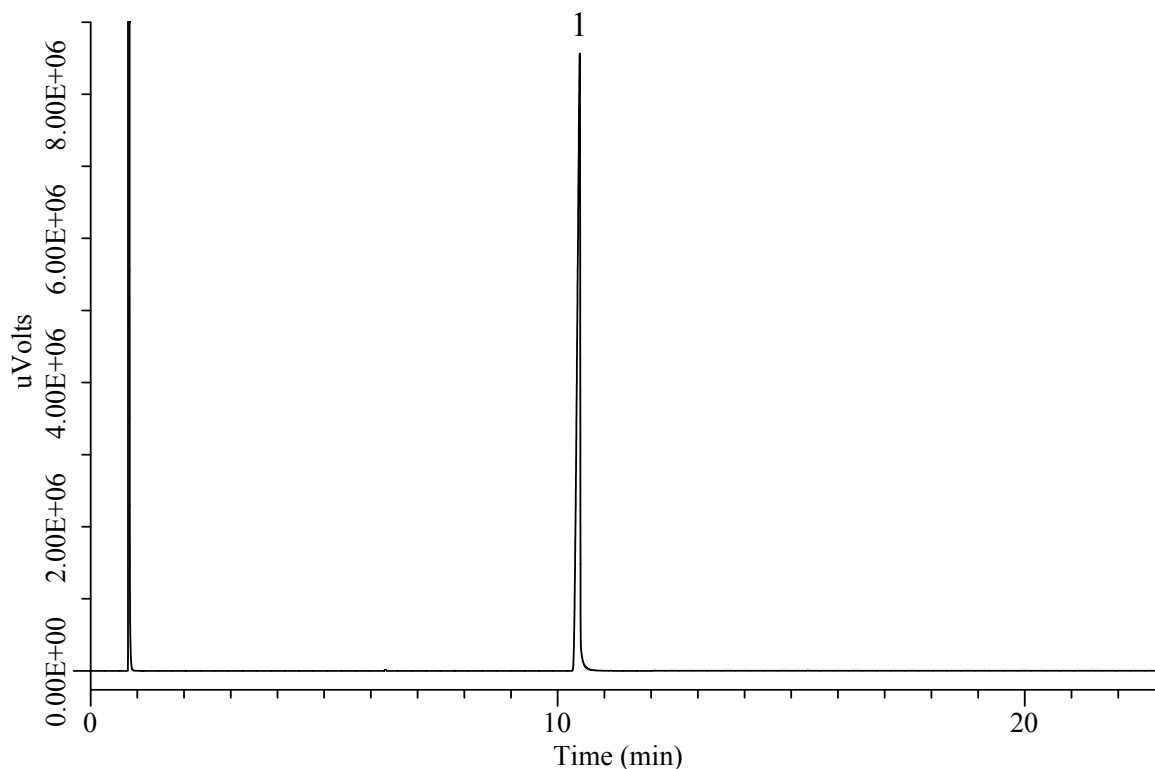

### Conditions

**System** : GC-FID  
**Column** : InertCap Pure-WAX  
0.53 mm I.D. x 30 m df = 1.0 µm  
**Col. Cat. No.** : 1010-68445  
**Col. Temp.** : 80 °C(1 min) - 10 °C/min - 190°C(20 min)  
**Carrier Gas** : He 13.0 mL/min Constant Flow  
**Injection** : Split 1:20  
200 °C  
**Injection Vol.** : 1 µL  
**Detection** : FID Auto Range  
200 °C  
**Sample** : Standard  
100 mg/mL in Methanol

### Analyte

1. 2-Pyrrolidone