


Aromatic Organic acids

Data No. LA106-0000

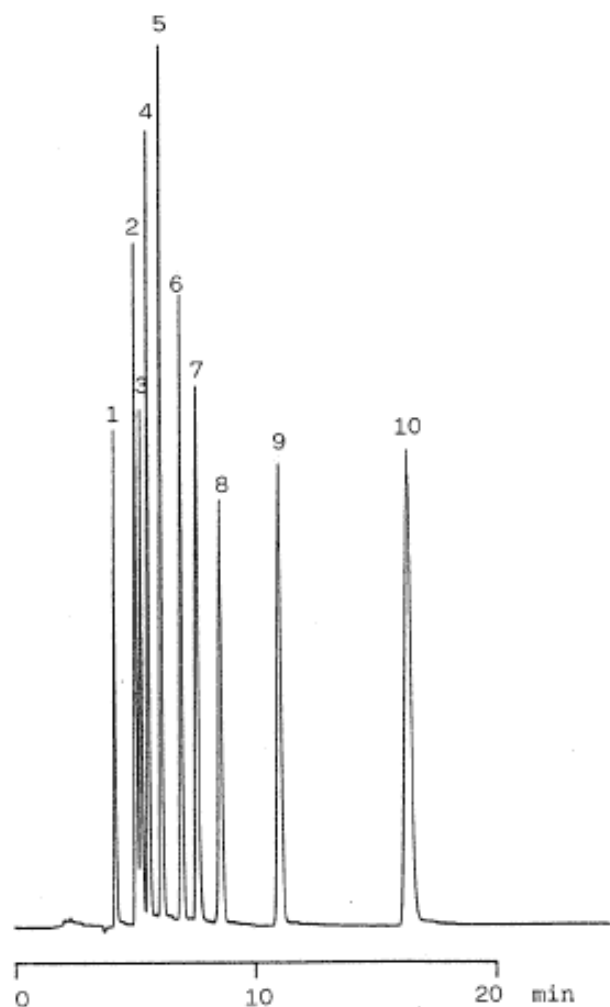

Column : Inertsil C8
5 μ m
250 \times 4.6 mmI.D.

Column Cat. No. : 5020-01202

Eluent : CH₃CN/0.1% H₃PO₄=25/75

Flow rate : 1.0 mL/min.

Col. Temp. : 40°C

Detection : UV 254nm

Injection Vol. :

Sample :

Solutes : 1. 3,4-Dihydroxy phenyl acetic acid
2. Caffeic acid
3. p-Hydroxyphenyl acetic acid
4. Vanillic acid
5. Gentisic acid
6. p-Coumaric acid
7. Ferulic acid
8. m-Coumaric acid
9. o-Coumaric acid
10. Salicylic acid