


# InertSearch™ for GC

InertCap® Applications

## Isomers of cresol

Data No. GA083-0000

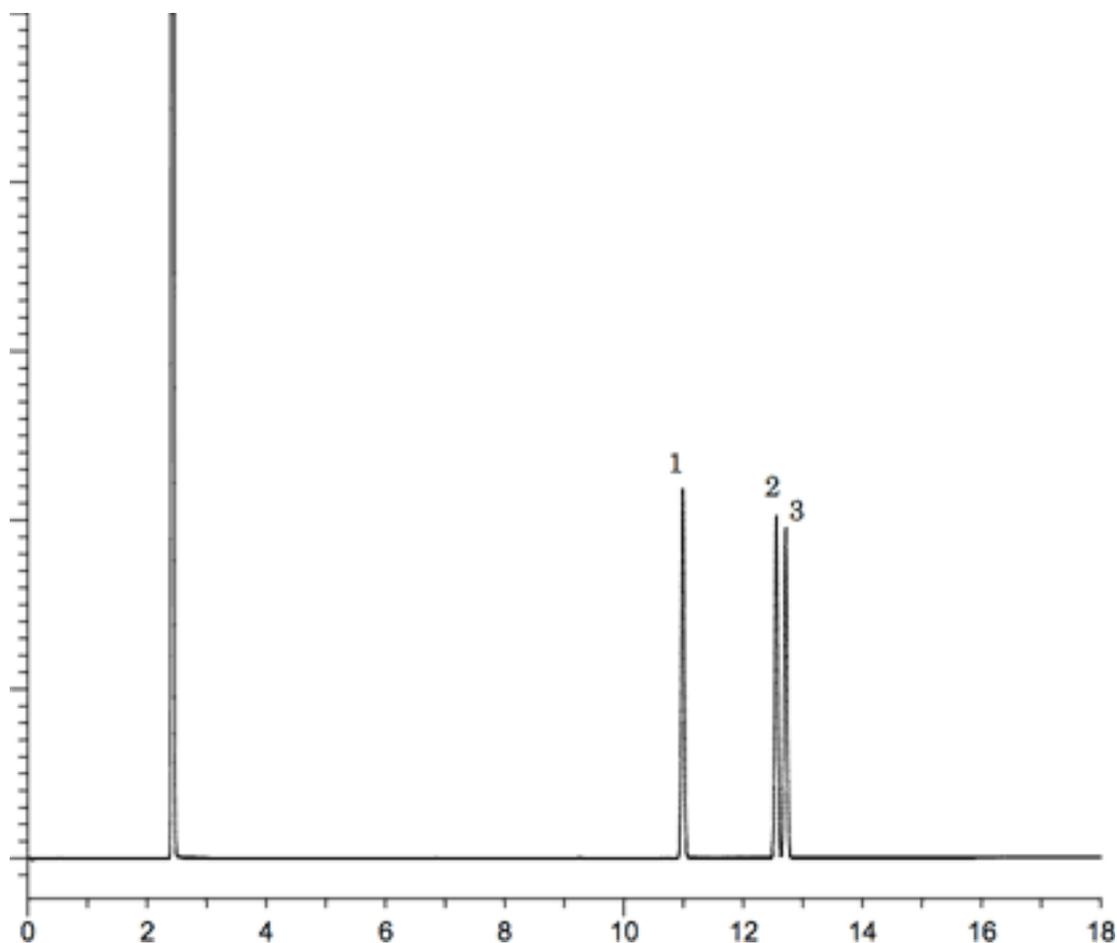

System : GC/FID  
Column : InertCap Pure-WAX  
0.25mm I.D. × 30m df=0.25µm  
Column Temp. : 150°C (5min) → 5°C/min → 250°C  
Carrier Gas : He 100kPa  
Injection : Split 1 : 50 0.2µl  
250°C  
Detection : FID Range 10°  
260°C  
Sample : Cresol isomers

Solutes : 1. o-Cresol  
2. p-Cresol  
3. m-Cresol