


# InertSearch™ for LC

Inertsil® Applications

## Analysis of Ergosterol(InertSustain C18)

Data No. LB014-0811

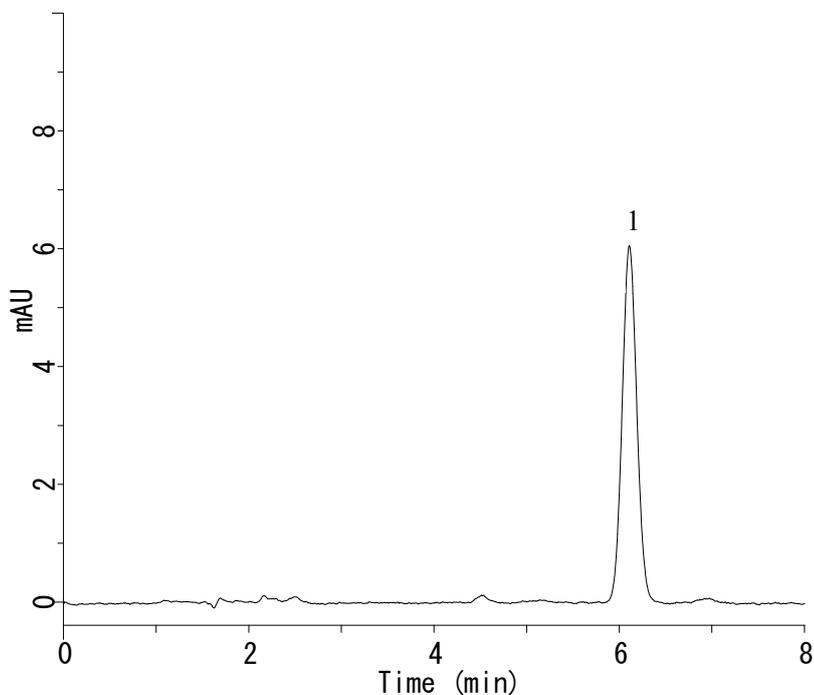

### Conditions

**System** : GL-7400 HPLC system  
**Column** : InertSustain C18 (5 µm, 150 x 4.6 mm I.D.)  
**Column Cat. No.** : 5020-07345  
**Eluent** : CH<sub>3</sub>OH  
**Flow Rate** : 1.0 mL/min  
**Col. Temp.** : 40 °C  
**Detection** : UV 290 nm (GL-7450 UV Detector)  
**Injection Vol.** : 10 µL  
**Sample** : Standard

### Analyte

1. Ergosterol (10 mg/L)