


# InertSearch™ for GC

InertCap™ Applications

## Flavor by Japanese Standards of Food Additives

Data No. GA316-0984

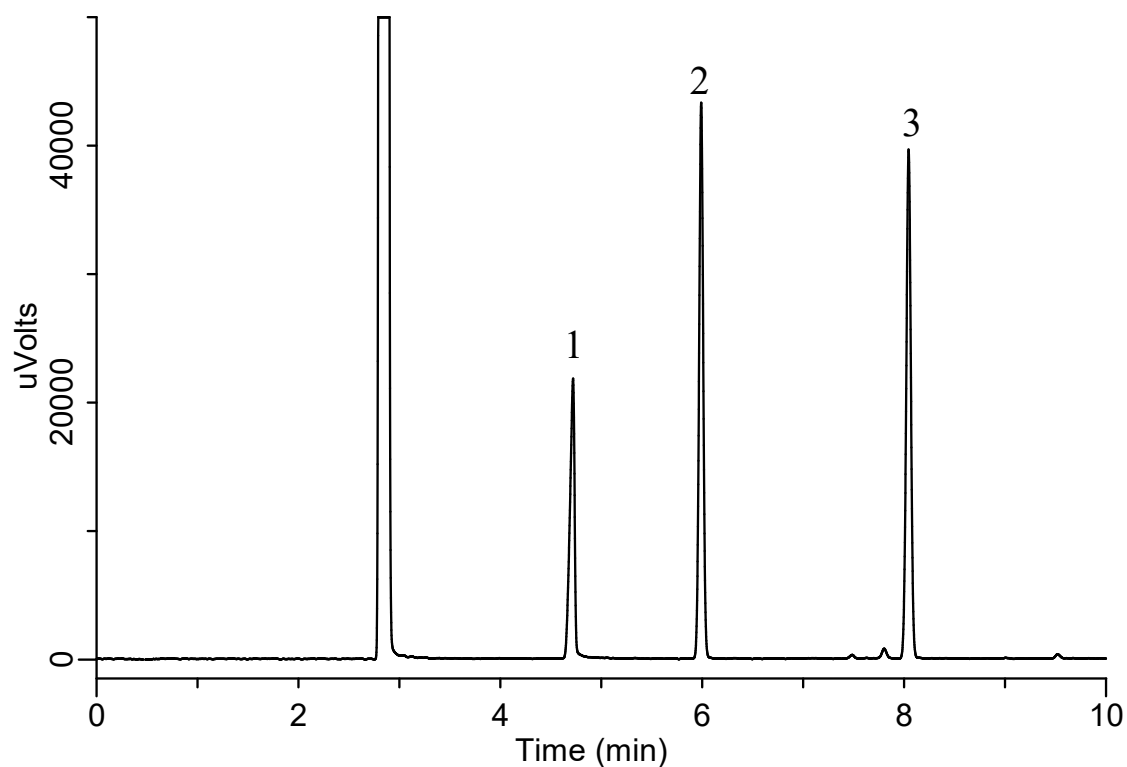

### Conditions

**System** : GC-FID  
**Column** : InertCap 1  
0.25 mm I.D. × 30 m df = 0.25 μm  
**Col. Cat. No.** : 1010-11142  
**Col. Temp.** : 100 °C - 5 °C/min - 230 °C (10 min hold)  
**Carrier Gas** : He 60 kPa  
**Injection** : Split flow 65 mL/min  
275 °C  
**Detection** : FID Auto Range  
300 °C  
**Sample Size** : 1000 μg/mL in Methanol  
1 μL

### Analyte

1. Hexanoic acid  
2. Acetophenone  
3. *p*-Methylacetophenone