

Endodontic Equipment








Woodpecker Apex Locator, with more than

accurate value to help you in workinglength determination during root canal treatment.

- ▲ 6th generation apex locator for measurement in presence of blood, pus, residual pulp or saliva etc.
- 97.71% accurate than the conventional radiography method
- Multi-frequency impedance technology precisely measure irrespective of the tooth type
- ★ File clip, hook, file probe can be easily Autoclavable to avoid crosscontamination
- ▲ Display: Multi-color 4.5" LCD display
- ▲ Main unit weight: 385gm
- ▲ Dimensions: 138mm x 77mm x 107mm



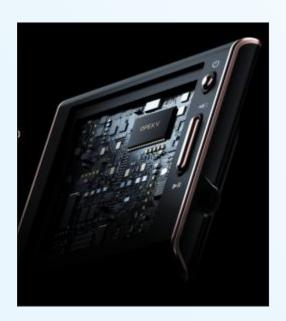

New Digital Circuit Industry recognized U.S. central processing chip

- ▲ Information processing speed increased by
- ▲ Adopt high-precision components imported from Japan
- ▲ High resolution signal without distortion
- ▲ Automatic calibration during 0.5s booting





Wave Q

Sonic Irrigation Activator

Wave Q is a specially designed Endo Activator for dentists. With the working tip, the flushing solution is activated through the working tip's vibration to clean the root canal.

13,000

vibrations per minute

62g

holding weight

90 days

battery duration

20mm

holding size











SIZE	176 X 22 X 22mm
WEIGHT	60±2G
BATTERY	3.7V / 800mAh lithium battery
OUTPUT	5V DC, 2A
INPUT	AC 100-240~50/60HZ
MATERIALS	Contra angle: plastic and silicone Working tip: POM plastic

Features

Easy to use

3 different sonic cleaning frequencies. Easy to install and disassemble in one step

Avoid breakage

Polymer material (POM) material working tips avoid breakage during cleaning.


Ergonomic design

Only 62g and the grip part diameter is 20mm It is the most comfortable size for holding.

Safe and hygienic

Steam sterilizable contra angle avoids cross-infection Sterilization condition: 135°C, 4 minutes, 0.22Mpa



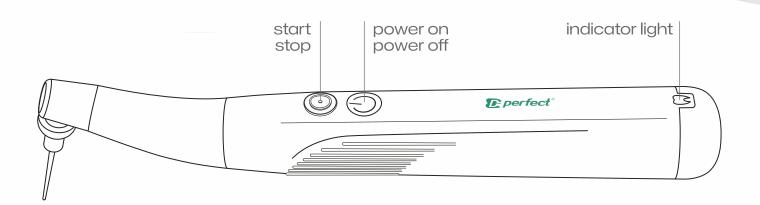

For any Sales & Support



Corp. Off.: 3 Local Shopping Centre, MOR Land, Near J Block, DDA Market, New Rajender Nagar, New Delhi -60





STREAM

Traditional anesthetic procedure is manual controlled and the extrusion speed is unstable, which made the poor absorption effect and the painful feeling of patient.

Stream X is an assistor booster which is specially designed in length and weight, the length is small and the center of gravity of the product is in the front. It relieves patients' pain and doctors' fatigue.

MODE	TIME
Fast Mode	1min
Slow Mode	2min
PDL Mode	4min
Custom Mode	1-4min



ASSISTOR BOOSTER





STREAM>

FEATURES

Custom Injection Setting

Stream X is equipped with 4 modes of injection time which varys from 1 to 4 minutes.

Short Reflow Setting

The reflow time is just 30 seconds. It's really convenient for the doctor to start the next operation.

Comfortable to Hold

Ergonomic & comfortable designed, weight 72g, holding size 20mm, the length is small and the center of gravity is in the front.

Durable and Efficient

The 350mAh lithium battery can be continuous used for 30 days (6 root canals per day).

SPECIFICATION

SIZE 166x21x33mm

WEIGHT 73±2g


BATTERY 3.7V/350mAh lithium battery

BATTERY DURATION 30 days

OUTPUT POWER 5V DC. 2A

INPUT POWER 100-240~50/60Hz

2h(Recommended to charge for more than 3 hours for the first time.) **CHARGING TIME**

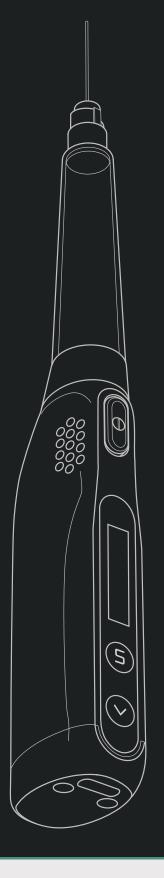

For any Sales & Support



Corp. Off.: 3 Local Shopping Centre, MOR Land, Near J Block, DDA Market, New Rajender Nagar, New Delhi -60







Gutta Percha Cutting Pen







GP FLOW

Gutta Percha Cutting Pen

<u>Features</u>

Efficient

Rapid heating in 0.2s at 4 levels temperature: 150°C, 180°C, 200°C, 230°C

Durable

Equipped with 2600mAh large battery and can be continuous used for 30 days

Easy to use

The three-button operation logic is simple, making it easy to use.

Comfortable

Ergonomic design: comfortable holding weight 95g and size 20mm, and the working tip is 360° rotatable.

Information

SIZE 174 x 25 x 25mm

WORKING TIP SIZE 40/.04; 45/.04; 50/.04; 55/.08

HEATING TEMPERATURE 150°C, 180°C, 200°C, 230°C

WEIGHT 95±2a

BATTERY 3.7V / 2600mAh lithium battery

OUTPUT POWER SUPPLY 5V DC, 2A


INPUT POWER SUPPLY AC 100-240~50/60Hz

STERILIZATION CONDITIONS 135°C, 4minutes, 0.22Mpa



For any Sales & Support



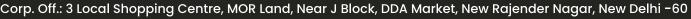





WHERE TRUST COMES FIRST



Warranty Ke Saath Bhi Warranty Ke Baad Bhi

For any Sales & Support

- 9599281852, 8448398852
- ≥ sales@unicorndenmart.com

Customer Care

011-45551250, 95994 48518

care@unicorndenmart.com | care1@unicorndenmart.com



SCAN TO VISIT OUR WEBSITE



<u>Corporate Office:</u> 3 Local Shopping Centre, MOR Land, Near J Block DDA Market, New Rajender Nagar, New Delhi – 110060







