


## Recycle preparative purification of Steroids

Data No. LA141-0000

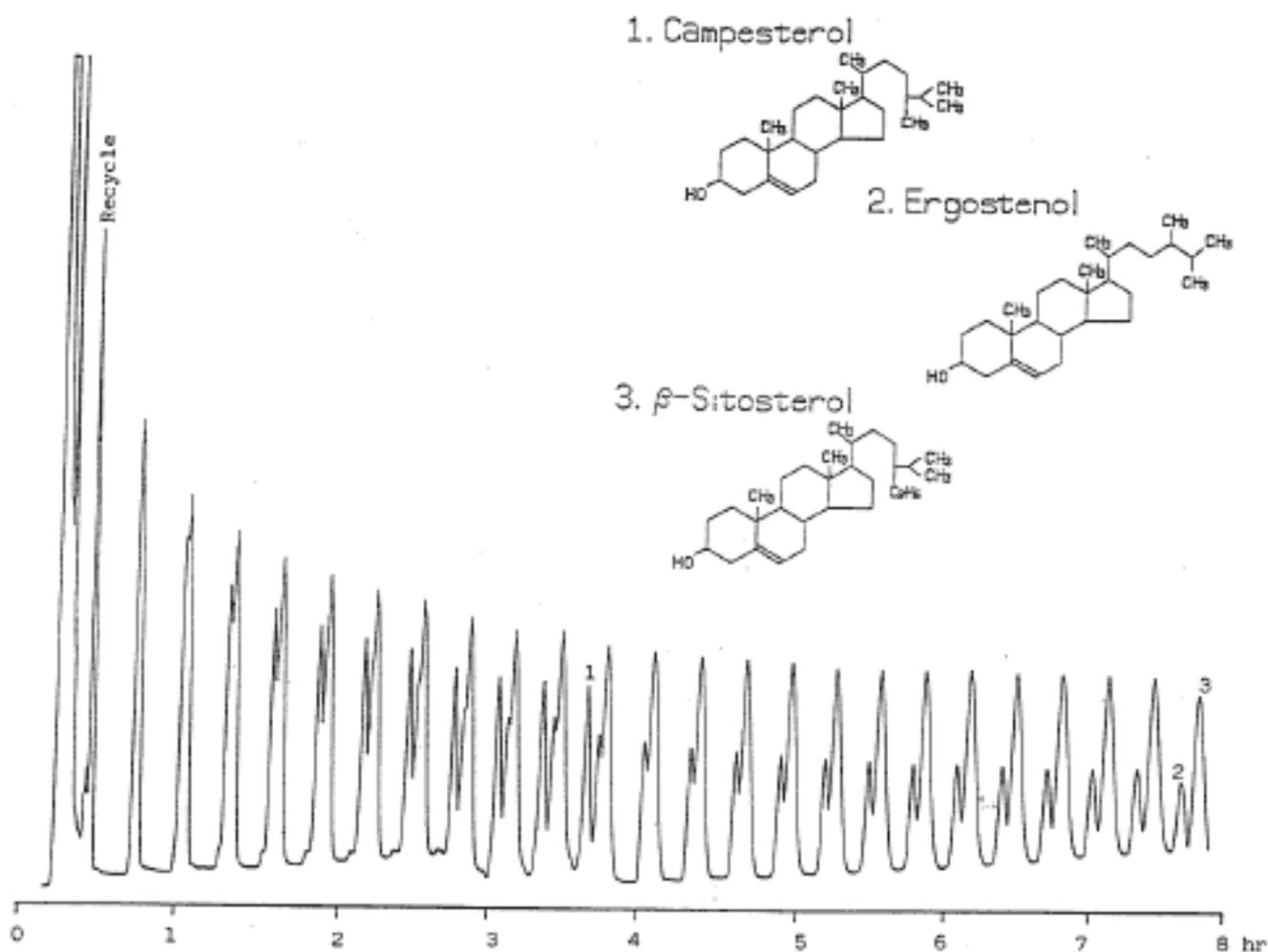

Column : Inertsil PREP-ODS  
10 $\mu$ m  
250  $\times$  20mm I.D.

Solutes : 1. Campesterol  
2. Ergosterol  
3.  $\beta$ -Sitosterol

Column Cat. No. : 5020-31123

Eluent : C<sub>2</sub>H<sub>5</sub>OH/H<sub>2</sub>O=93/7

Flow rate : 12 mL/min.

Col. Temp. : 40°C

Detection : UV 210 nm

Injection Vol. : 2mL

Sample : Crude  $\beta$ -sitosterol(20mg/mL)