


Analysis of Xylene Isomers

Data No. LA477-0000

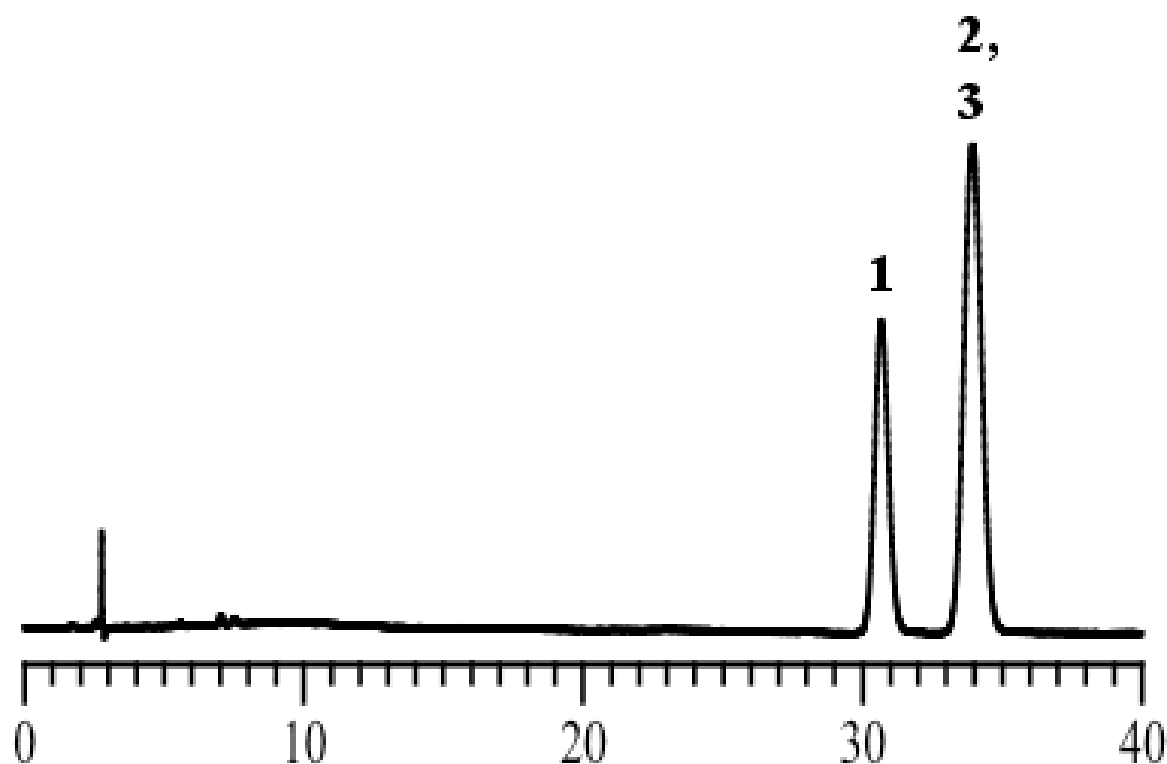

Column : Inertsil ODS-3
5 μ m
250 x 4.6 mm I.D.

Solutes : 1. o-Xylene
2. m-Xylene
3. p-Xylene

Column Cat. No. : 5020-01732

Eluent : A) Ethanol
B) β -Cyclodextrin/H₂O = 5/995, v/v
A/B = 50/50, v/v

Flow rate : 1.0 mL/min.

Col. Temp. : 40°C

Detection : UV 254 nm

Injection Vol. :

Sample :