


Analysis of Catecholamines

Data No. LB010-0721

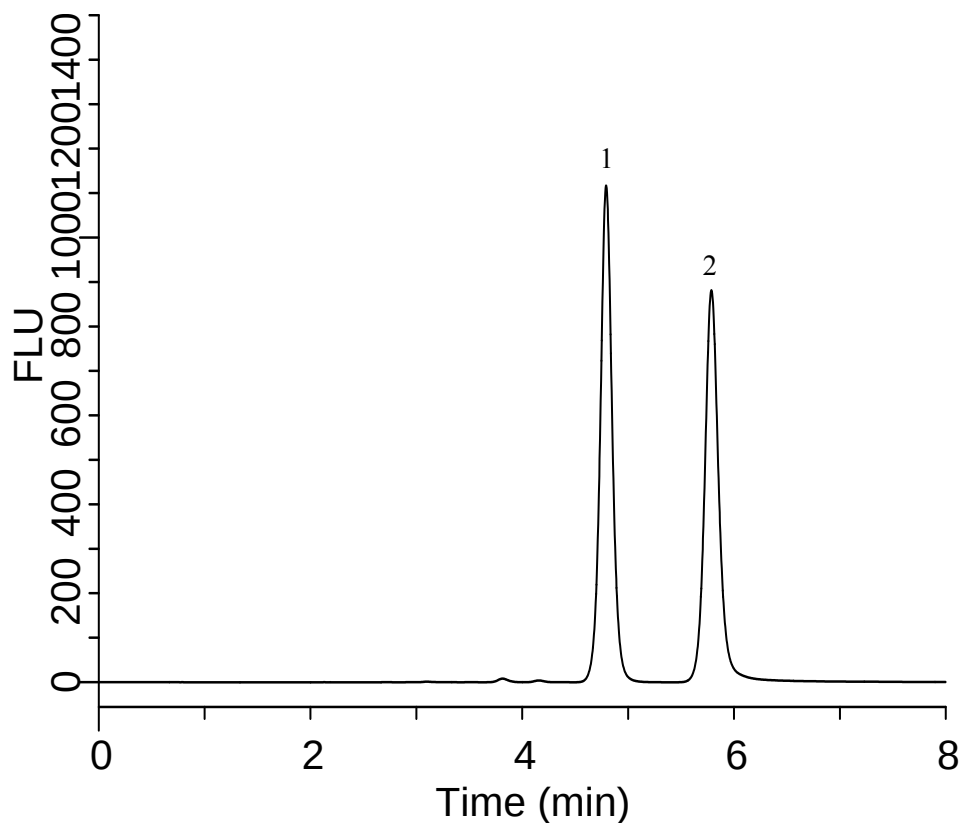

Conditions

System : GL-7400 HPLC system
Column : Inertsil Amide (5 μ m, 250 x 4.6 mm I.D.)
Column Cat. No. : 5020-07836
Eluent : A) CH₃CN
B) 20 mM HCOONH₄ (pH 3.0, HCOOH)
A/B = 65/35, v/v
Flow rate : 1 mL/min
Col. Temp. : 40 °C
Detection : FL Ex 280 nm Em 350 nm (GL-7453A FL Detector)
Injection Vol. : 10 μ L
Sample : Standard

Analyte:

1.	Normetanephrine	100 μ g/L
2.	Norepinephrine	100 μ g/L