


InertSearch™ for LC

Inertsil® Applications

Analysis of Gamma-Oryzanol

Data No. LA249-0000

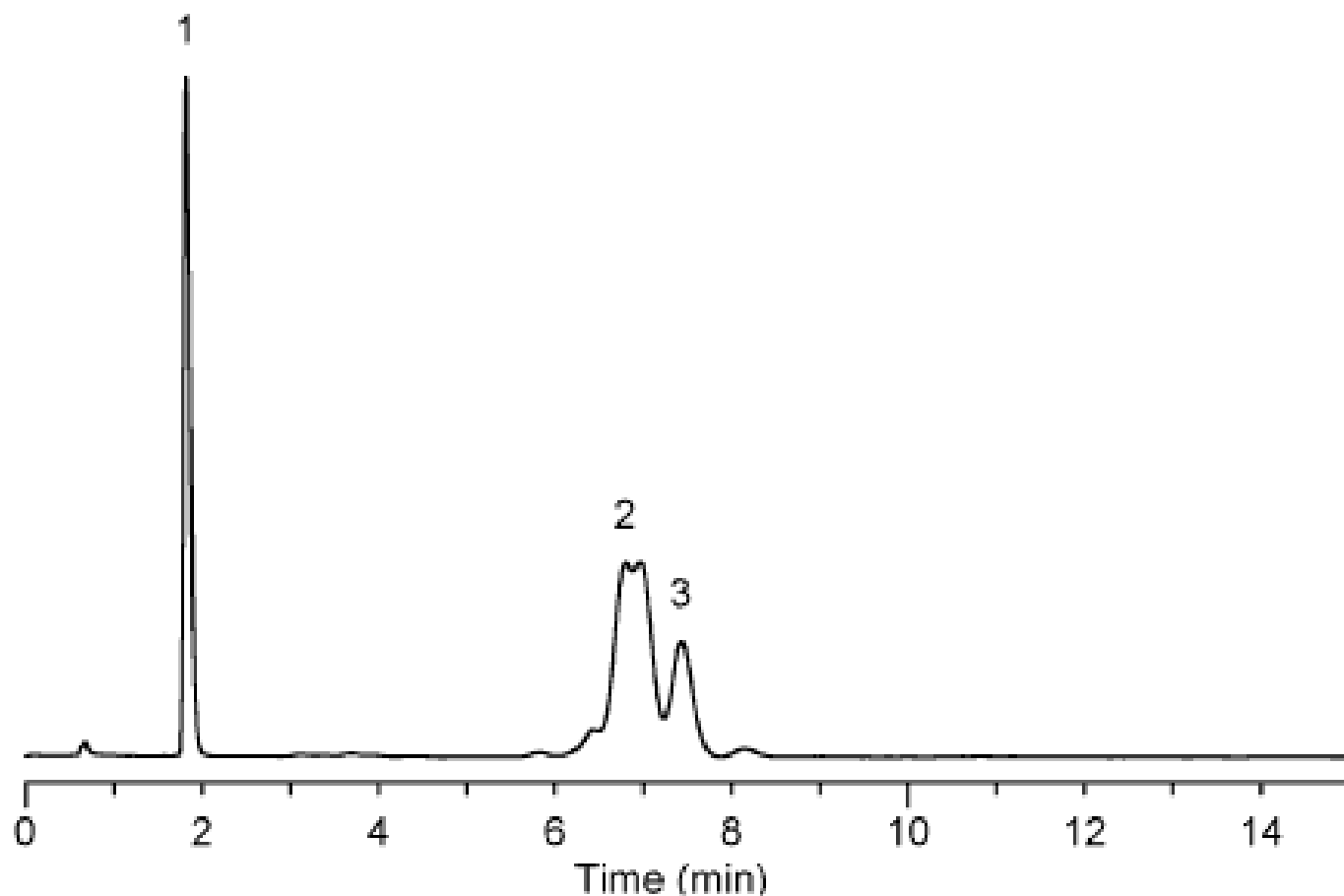

Column	: Inertsil C8-3	Solutes	: 1. vanillin (50ug/mL)
	5µm		2. γ -oryzanol (mixture 100ug/mL)
	150 × 4.6mm I.D.		
Column Cat. No.	: 5020-01900		
Eluent	: Methanol: Water=95:5		
Flow rate	: 1.0mL/min		
Col. Temp.	: 40°C		
Detection	: UV 320 nm		
Injection Vol.	: 1µL		
Sample	: γ -Oryzanol		