


InertSearch for LC

Inertsil Applications

Analysis of Glutathione

Data No. LB730-0806

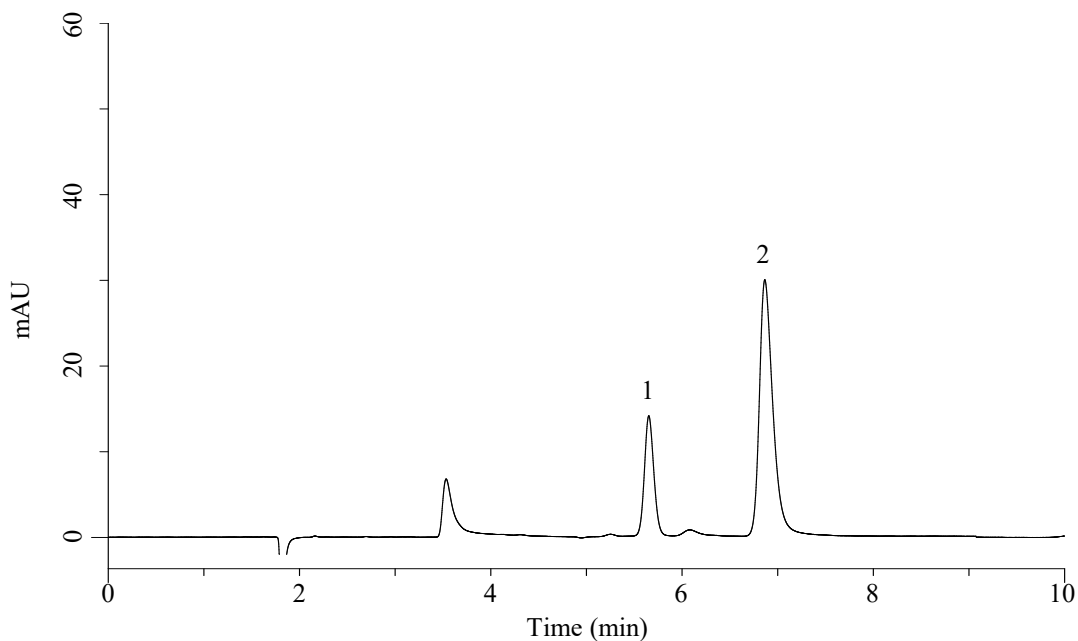

Conditions

System : Chromaster HPLC system (HITACHI)
Column : InertSustain AX-C18 (GL Sciences Inc.)
(3 μ m, 150 x 2.1 mm I.D.)
Column Cat. No. : 5020-91039
Eluent : 10 mM HCOONH₄ in H₂O
Flow Rate : 0.2 mL/min
Col. Temp. : 40 °C
Detection : UV 220 nm
Injection Vol. : 1 μ L
Sample : Standard

Analyte:

1. Glutathione (reduced form)
2. Glutathione (oxidized form)
(500 mg/L each)