


Analysis of Gamma-Oryzanol

Data No. LA248-0000

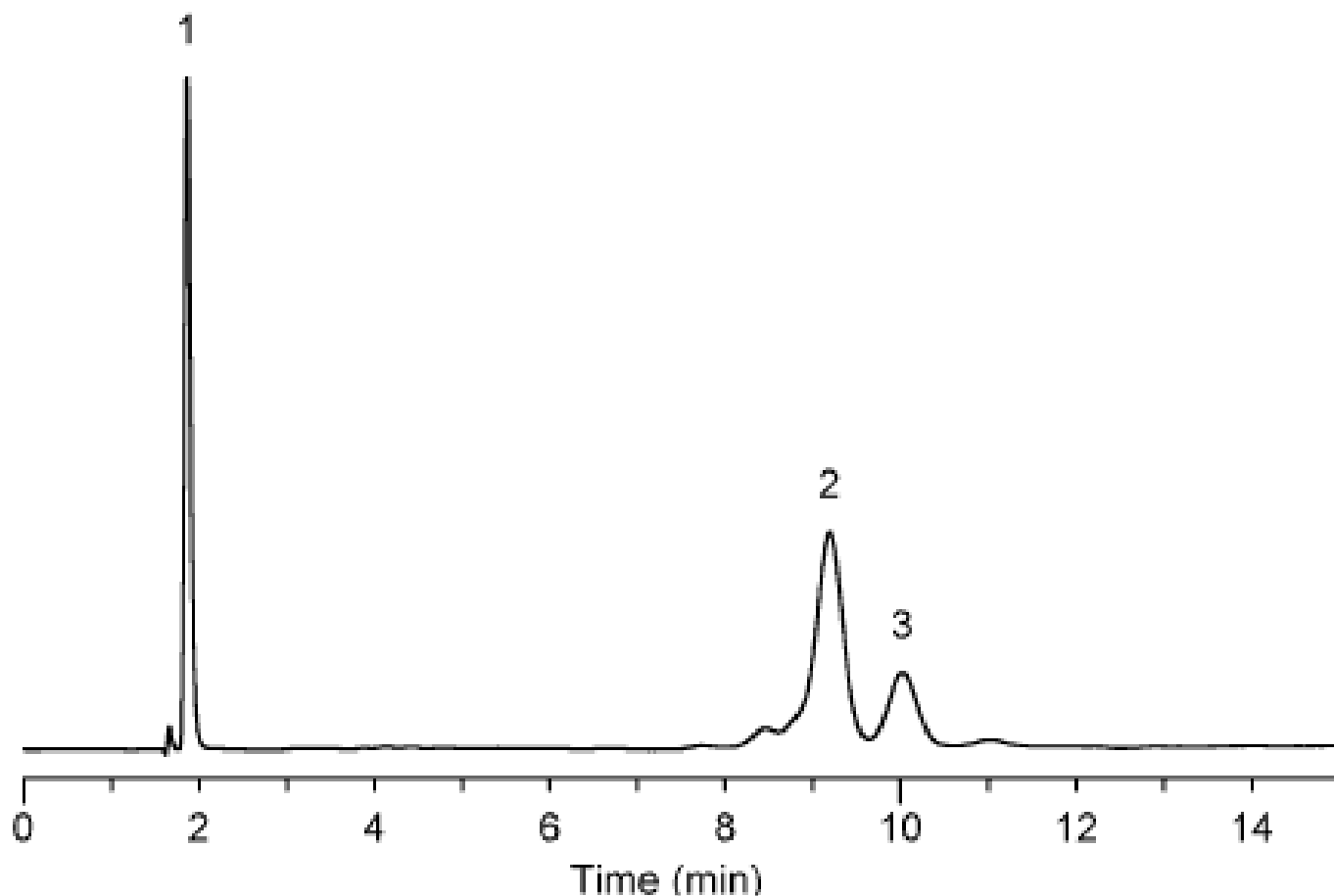

Column	: Inertsil ODS-3V	Solutes	: 1. vanillin
	5 μ m		2. γ -oryzanol-1
	150 \times 4.6mm I.D.		3. γ -oryzanol-2

Column Cat. No. :

Eluent : Methanol:2-Propanol:Water=70:30:1

Flow rate : 1.0mL/min

Col. Temp. : 40°C

Detection : UV 320 nm

Injection Vol. : 1 μ L

Sample : γ -Oryzanol