


InertSearch™ for GC

InertCap® Applications

Organic solvents

Data No. GA097-0000

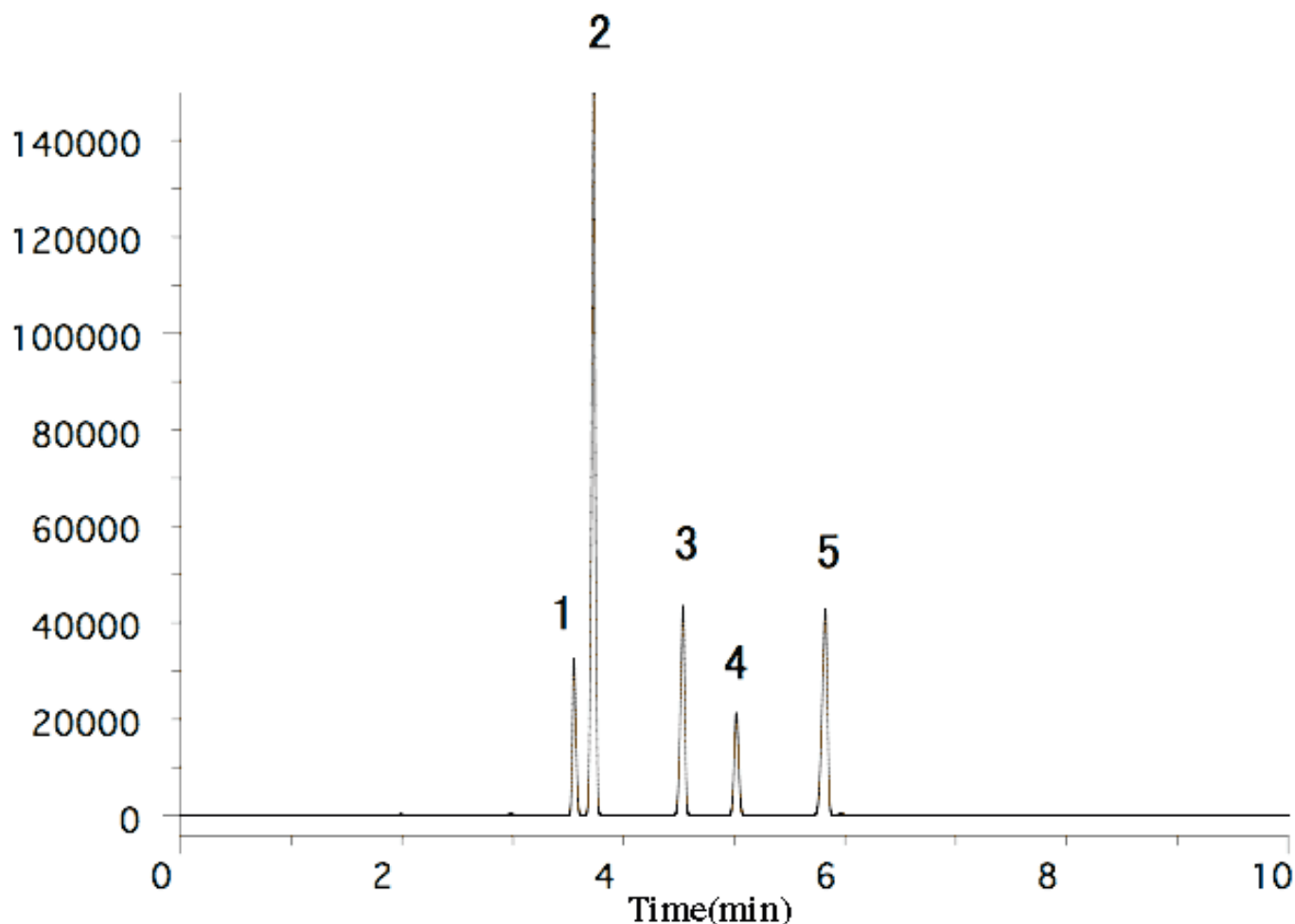

System : GC/FID
Column : InertCap Pure-WAX
0.53mm I.D. × 30m df=1µm
Column Temp. : 40°C – 5°C/min – 90°C
Carrier Gas : He 20kPa 35cm/sec
Injection : Split Flow 200ml/min
240°C
Detection : FID Range10²
240°C
Sample : Solvents
Mixed evenly
Sample Size : 0.4µL

Solutes : 1. Methylene chloride
2. Benzene
3. Trichloroethylene
4. Chloroform
5. 1,4-Dioxene