


## Analysis of Nucleic Acid Bases

Data No. LA250-0000

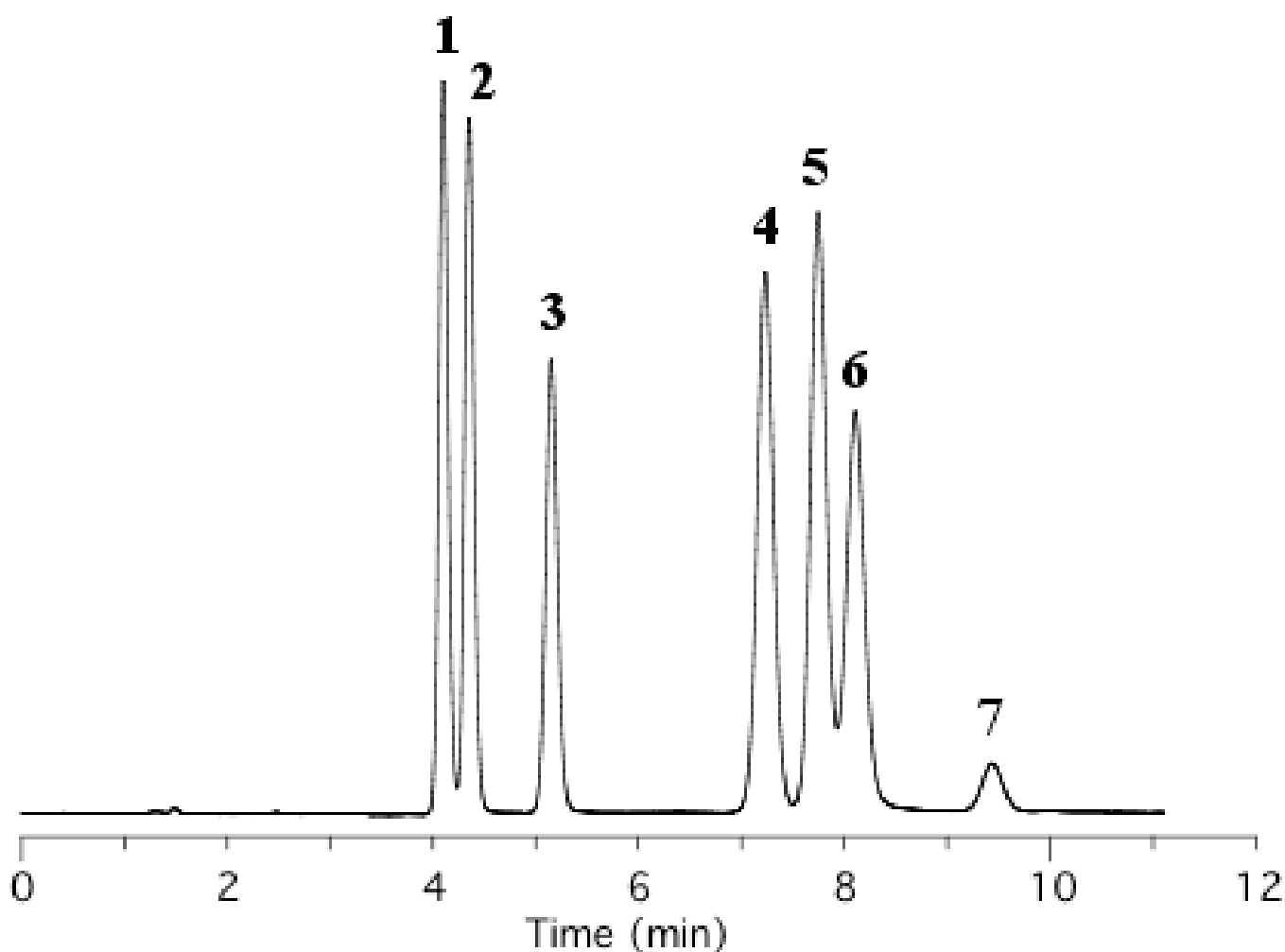

Column	: Inertsil Diol 5 µm 250 x 4.6 mm I.D.	Solutes	: 1. Thymine 2. Uracil 3. Uridine 4. Adenosine 5. Cytosine 6. Cytidine 7. Guanosine
Column Cat. No.	: 5020-05646		
Eluent	: A) CH <sub>3</sub> CN B) H <sub>2</sub> O A/B = 85/15, w/w		
Flow rate	: 1.0 mL/min.		
Col. Temp.	: 40°C		
Detection	: UV 254 nm		
Injection Vol.	:		
Sample	:		