


InertSearch™ for GC

InertCap® Applications

Styrene, Xylenes

Data No. GA082-0000

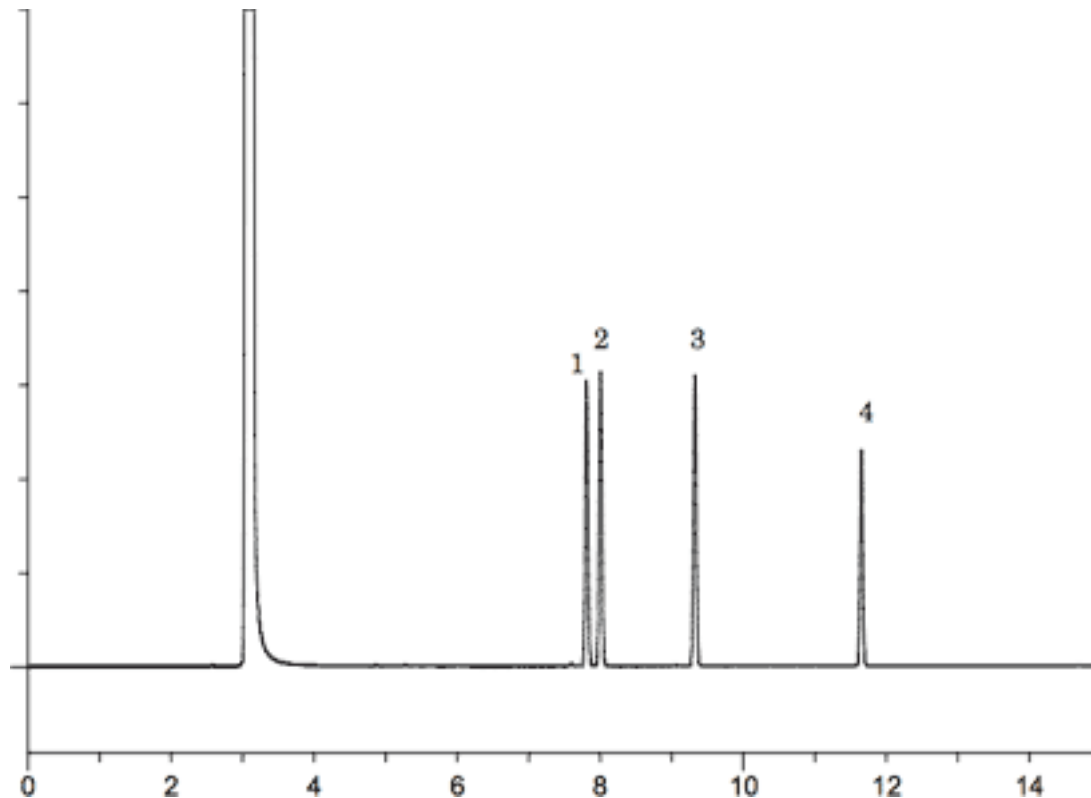

System : GC/FID
Column : InertCap Pure-WAX
0.25mm I.D. x 30m df=0.25µm
Column Temp. : 50°C (5min) → 5°C/min → 200°C
Carrier Gas : He 100kPa
Injection : Split 1 : 50 0.2µl
250°C
Detection : FID Range 10°
260°C
Sample : Styrene, Xylenes
1000mg/l Solvent : Methanol

Solutes : 1. p-Xylene
2. m-Xylene
3. o-Xylene
4. Styrene