


Phthalic acid esters

Data No. GA176-0904

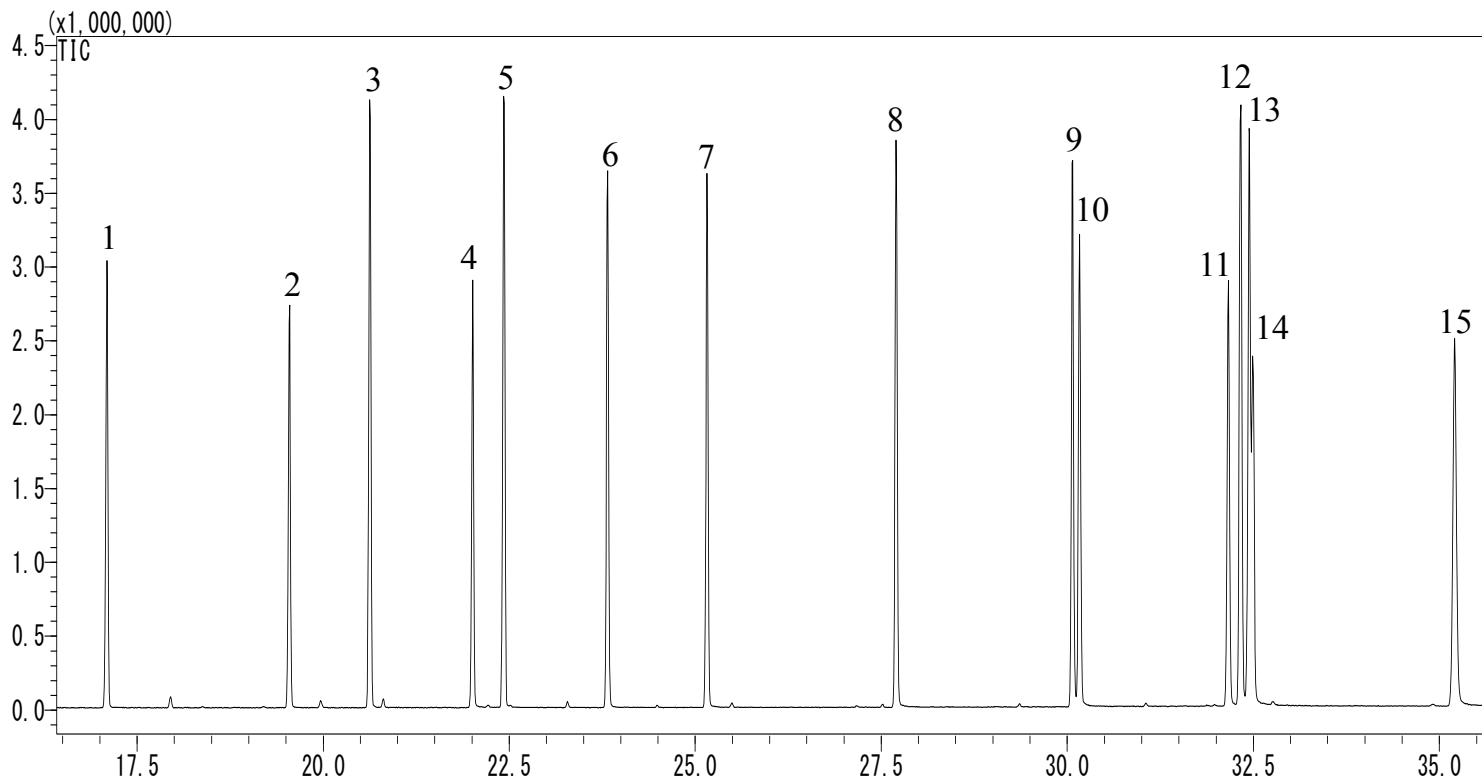

Conditions

System : GC/MS
Column : InertCap 5MS/NP
0.25 mm I.D. x 30 m df = 0.25 µm
Col.Cat. No. : 1010-18642
Col. Temp. : 60 °C(3 min hold) - 7 °C/min - 260 °C(19 min hold)
Carrier Gas : He 36.5 cm/sec
Injection : Splitless 1 min
250 °C
Detection : MS Scan (35 – 500 m/z)
Interface Temp. 280 °C
Sample Size : 10 µg/mL in Acetone
1 µL

Analyte :

1. Dimethyl phthalate
2. Diethyl phthalate
3. Di-isopropyl phthalate
4. Diallyl phthalate
5. Di-*n*-propyl phthalate
6. Di-isobutyl phthalate
7. Di-*n*-butyl phthalate
8. Di-*n*-pentyl phthalate
9. Di-*n*-hexyl phthalate
10. Butylbenzyl phthalate
11. Dicyclohexyl phthalate
12. Di-*n*-heptyl phthalate
13. Di-2-ethylhexyl phthalate
14. Diphenyl phthalate
15. Di-*n*-octyl phthalate