


InertSearch for GC

InertCap Applications

Residual solvents in pharmaceuticals (JP18 Standard stock solution Class 2 C)

Data No. GA 333-0155

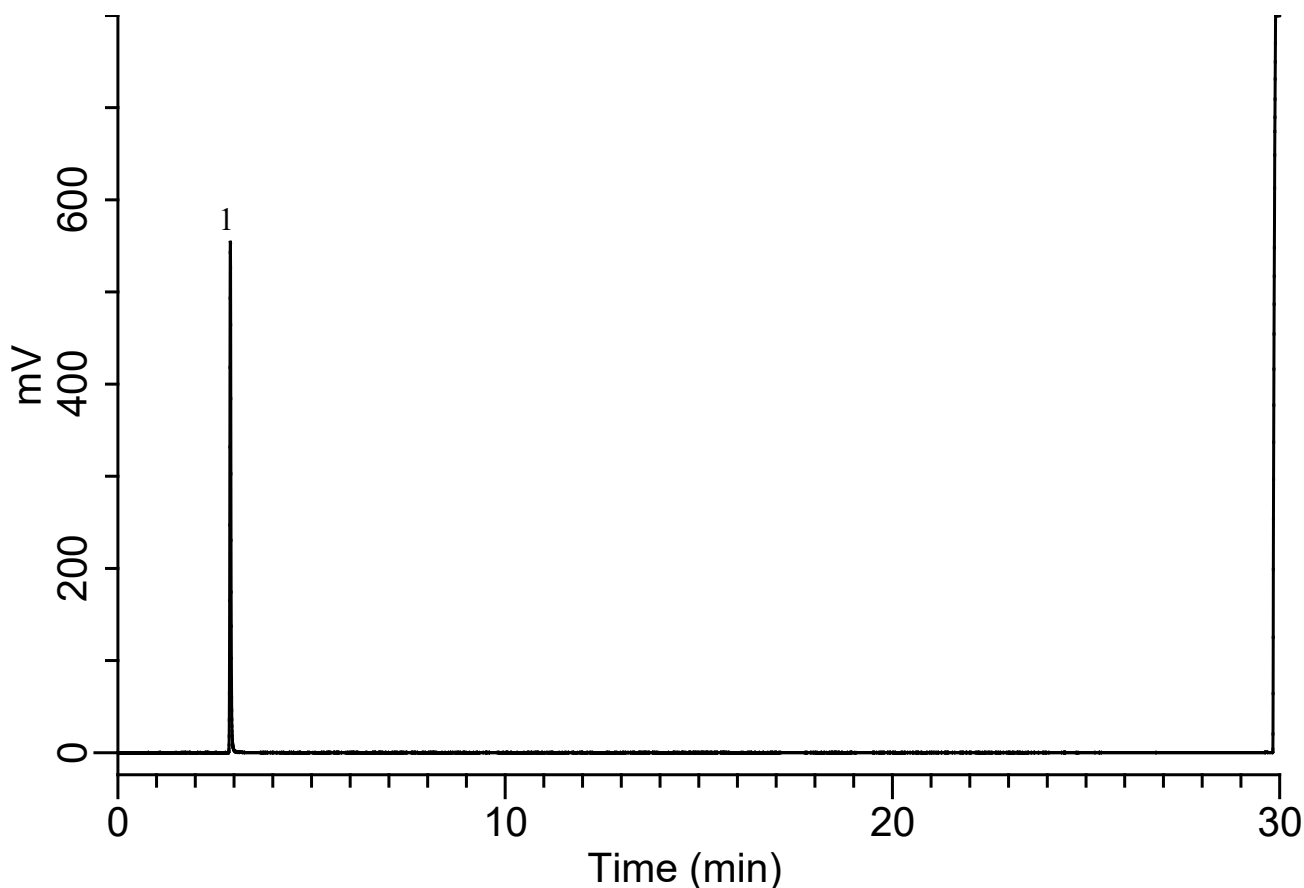

Conditions

System : GC4000 Plus-FID
Column : InertCap WAX
0.32 mm I.D. × 30 m df = 0.25 μm
Col. Cat. No. : 1010-67242
Col. Temp. : 50 °C(20 min) – 6 °C/min
- 165 °C(20 min)
Carrier Gas : He 2.20 mL/min (35cm/sec)
Injection : Split 20:1
250 °C
Detection : FID Auto Range
250 °C
Sample Size : 1.0 μL
Analyte in Dimethyl sulfoxide

Analyte (R.T. min)

1. Methyl isobutyl ketone (2.903)