


# InertSearch™ for GC

InertCap® Applications

## Sterols

Data No. GA032-0000

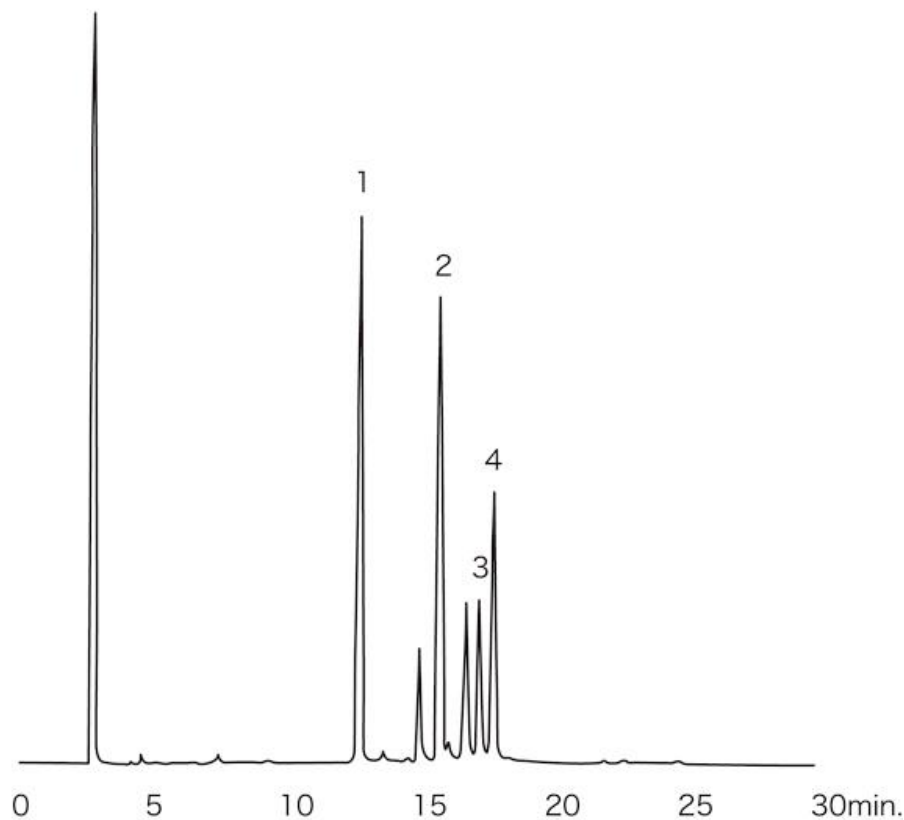

System : GC/FID  
Column : InertCap 1  
0.25mm I.D.×30m df=0.4µm  
Column Temp. : 300° C  
Carrier Gas : He 260kPa  
Injection : Split 1 : 10 1µL  
Detection : FID Range 1×128  
Sample : Sterols

Solutes : 1. Cholesterol  
2. Stigmasterol  
3.  $\beta$ -Sitosterol  
4. Lanosterol