


InertSearch™ for GC

InertCap® Applications

Organophosphorus pesticides

Data No. GA009-0000

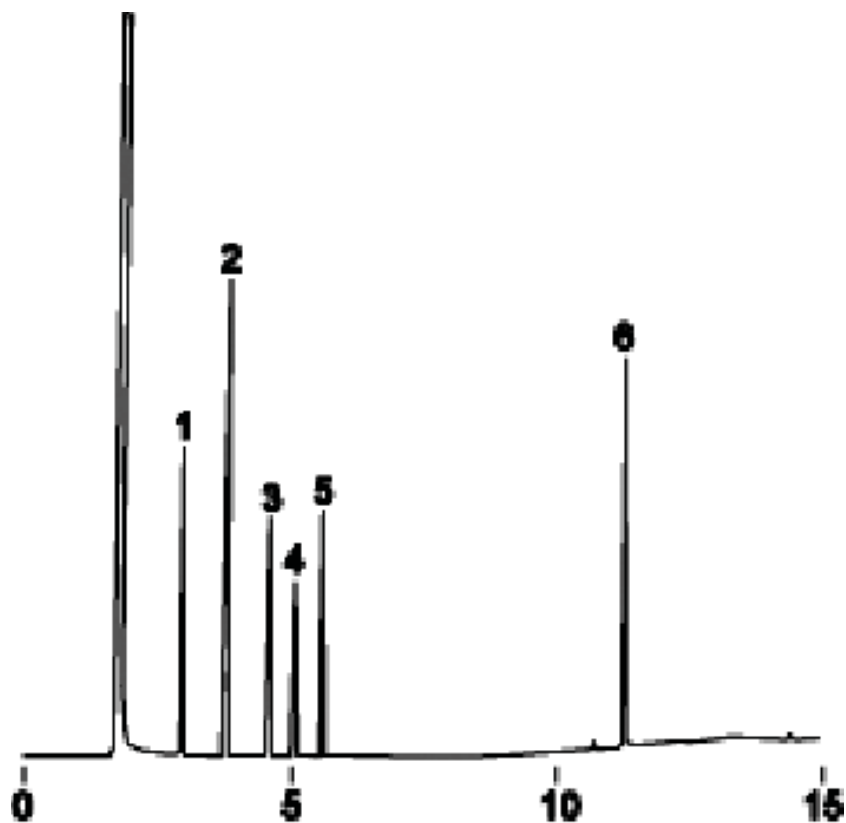

System : GC/FID
Column : InertCap 1
0.25mmI.D.×30m df=0.4µm
Column Temp. : 220° C (6 min)→10° C/min→280° C
Carrier Gas : He 80kPa
Injection : Split 1 : 50 1µl
Detection : FID Range 1×128
Sample : Pesticides

Solutes : 1. Methyl demetone
2. Diazinon
3. Methyl parathion
4. Sumithion
5. Parathion
6. EPN