


InertSearch™ for LC

Inertsil® Applications

Tetracyclines in Control Serum E28

Data No. LA155-0000

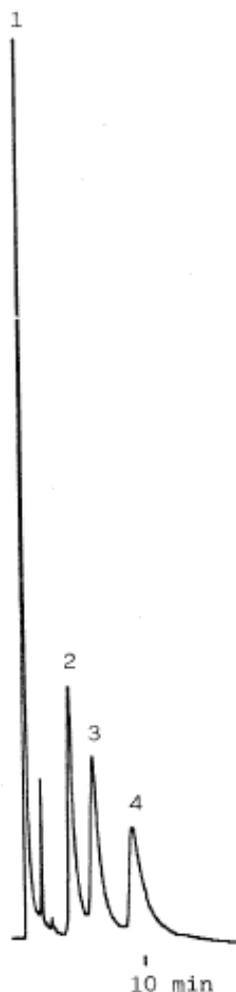

Column : Bioptic AV-2
150 × 4.6mm I.D.

Column Cat. No. : 5020-02512

Eluent : 50mM K₂HPO₄(pH7.5)/CH₃CN=95/5

Flow rate : 1.0 mL/min.

Col. Temp. : 40°C

Detection : UV 280 nm

Injection Vol. :

Sample :

Solutes : 1. Serum
2. Tetracycline
3. Doxycycline
4. Minocycline